

**Submittal Data Sheet**

**Features and Benefits**

The Frontall™ DISS medical gas recessed outlets represent the most technically advanced means of providing medical gas connections to the central piping system. The outlets incorporate a unique patented **front-loaded** feature which makes the Tri-Tech Medical outlets **easy to repair**, saving valuable time and money for the health care facility.

- **Five year parts and one year labor limited warranty\***
- UL Listed to UL 1331 and CSA Z9170-1.
- Conforms to NFPA 99 and CGA standards
- Accepts DISS 1220 gas specific nut & nipple adapter for Dental Vacuum and DISS 2080 gas specific nut and nipple adapter for Dental Air.
- Ease of maintenance without removing the faceplate using a cartridge insertion tool.
- Compensates for variation in wall thicknesses from 1/2" to 1 1/4"
- Modular design
- All outlets are 100% tested for flow and leaks
- 5" centers
- 360° swivel inlet tube
- Metal back and nameplate
- Strength of connections provides support for dispensing equipment without the use of additional brackets
- Made in the U.S.A.

\* See terms and Conditions, Document No. 99-0477, on our website at: [www.tri-techmedical.com](http://www.tri-techmedical.com) for complete details.

**Specification**

The medical gas wall outlets shall be Tri-Tech Medical Frontall™ DISS medical gas outlets. The outlets shall be modular singles, and shall incorporate nameplate color-coding. The single modular outlets shall be of a design that provides for ganging of rough-in plates in the field to form multiples. The gas services shall be sequentially arranged and located as shown on the plans with a minimum centerline spacing of 5 inches (12.5cm) between outlets.

The medical gas outlets shall be designed so that, once installed, routine service of both the primary and secondary check valves can be accomplished **without** removing the nameplate or gas specific portions. The primary check valve shall be unitized and of the cartridge type. So that the primary check valve can be removed for service without shutting off the gas supply to the outlet, the secondary check valve shall operate automatically to stop the free flow of the positive pressure gas. There shall be

no secondary check for vacuum or WAGD services. Medical gas outlets which require the removal of the nameplate or gas specific components for routine service shall not be



Shown above  
Dental Vacuum Part Number TDN91XF1122 and  
Dental Air Part Number TDN91XF112080


acceptable. The seal between the nameplate and the rough-in shall be a double O-ring which is serviceable without shutting off the gas service and is capable of supporting 10 pounds @ 2" from the nameplate without leaking.

The outlet nameplate shall be permanently color-coded with a durable scratch resistant and protective label. The outlet trim plate shall be durable plastic, attached with the nameplate to the rough-in assembly, and provide automatic plaster adjustment from 1/2 to 1 1/4 of an inch (1.3 cm to 3 cm). The name of the gas service shall be permanently marked on the outlet bracket and the chrome plated brass outlet body. The outlet's rough-in supply tube shall be a 7 inch (18 cm) length of 1/2 OD copper Type K, and labeled with the name of the gas service.

Medical gas outlets shall be cleaned for oxygen service in accordance with the current Compressed Gas Association (CGA) Pamphlet G-4.1, capped and placed in a protective container for shipment. The outlets shall be installed in strict accordance with manufacturer's instructions, and tested before use in conformance with local and federal codes.


MRI compatible models are available (See ordering information chart) and are made with non-ferrous materials.



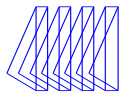

**Layout**

DIMENSION IN. (CM)

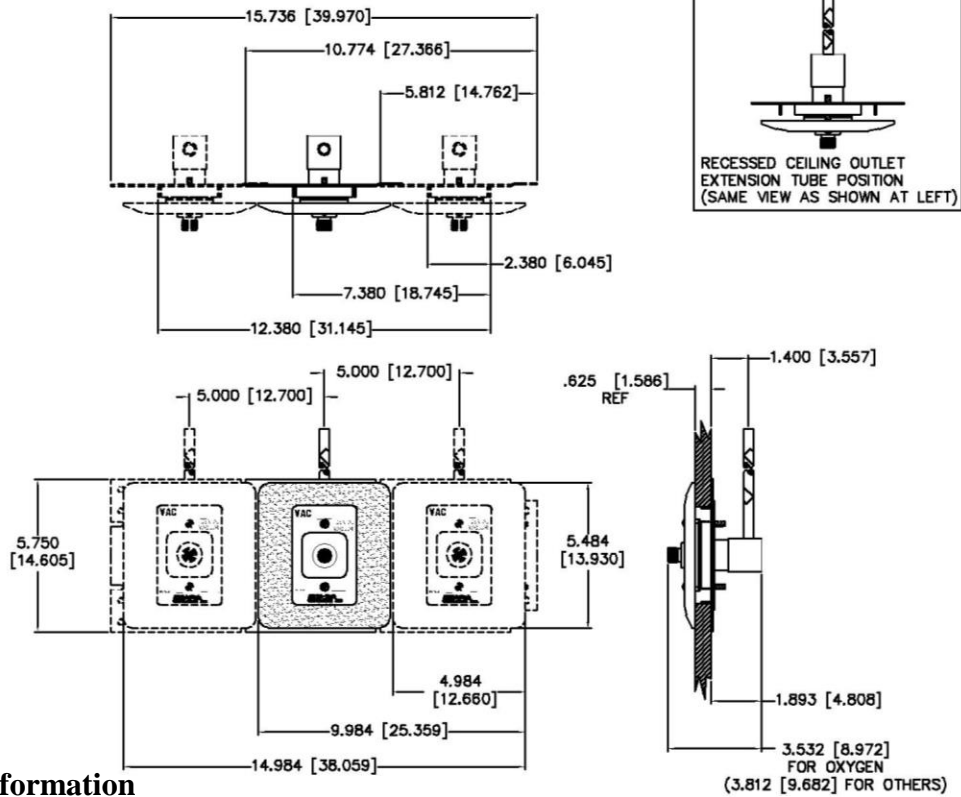

**Ordering Information**

*Standard Recessed Outlets Complete - with the Plastic White Trim Plate*

Gas Service	Recessed Complete		MRI Compatible	
	Wall	Ceiling	Wall	Ceiling
Dental Vacuum	TDN91XF1122	TDN91XF1222	TDN91XF1122-MRI	*TDN91XF1222-MRI
Dental Air	TDN91XF112080	TDN91XF122080	*TDN91XF112080-MRI	*TDN91XF122080-MRI
Slide	XF1144	N/A	XF1144-MRI	N/A

*\*Indicates that these are NOT UL listed Outlets  
Slides are not subject to UL Listing.*

*If Ordering Die Cast Trim Plates for Wall Outlets you need Back Bodies & Faceplates ordered Separately*

Gas Service	Back Body ½ Tube Size Wall	MRI Compatible Wall	Outlet Faceplate	MRI Compatible	Trim Plate Die Cast White	Trim Plate Die Cast Pewter
Dental Vacuum	TDN91XC1122-A	TDN91XC1122-A-MRI	TDN91XF1022-A	TDN91XF1022-A-MRI	XC1000-9DC	XC1000-9DCP
Dental Air	TDN91XC112080-A	*TDN91XC112080-A-MRI	TDN91XF102080-A	*TDN91XF102080-A-MRI		
Slide	Die Cast Pewter XF1144-DCP	Pewter MRI Compatible XF1144-DCP-MRI	Die Cast White XF1144-DC	Die Cast White MRI Compatible XF1144-DC-MRI		

*If Ordering Die Cast Trim Plates for Ceiling Outlets you need Back Bodies & Faceplates ordered Separately*

Gas Service	Back Body ½ Tube Size Ceiling	MRI Compatible Ceiling	Outlet Faceplate	MRI Compatible	Trim Plate Die Cast White	Trim Plate Die Cast Pewter
Dental Vacuum	TDN91XC1222-A	*TDN91XC1222-A-MRI	TDN91XF1022-A	TDN91XF1022-A-MRI	XC1000-9DC	XC1000-9DCP
Dental Air	TDN91XC122080-A	*TDN91XC122080-A-MRI	TDN91XF102080-A	*TDN91XF102080-A-MRI		

To place an order or request additional information call toll-free 1-800-253-8692 or fax us at 1-440-937-5060 or e mail us at [sales@tri-techmedical.com](mailto:sales@tri-techmedical.com). Please visit our web-site at [www.tri-techmedical.com](http://www.tri-techmedical.com)