


**Tri-Tech Medical Outlet**

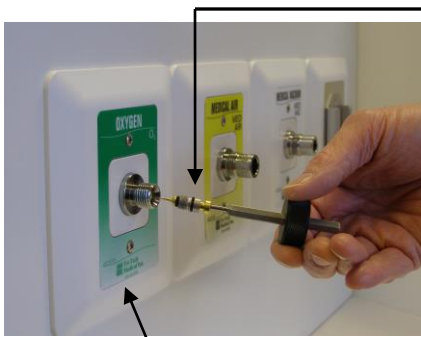


DISS Frontall Design – service outlet w/o removing nameplate (DISS only)

Durable, protective scratch resistant Vinyl label may be customized (i.e. – private labeling, color codes to meet various standards)

Brass Poppet machined surfaces  
Positive O-ring seal w/ mating surface

O-Rings can be serviced w/o shutting off- gas or removing secondary check

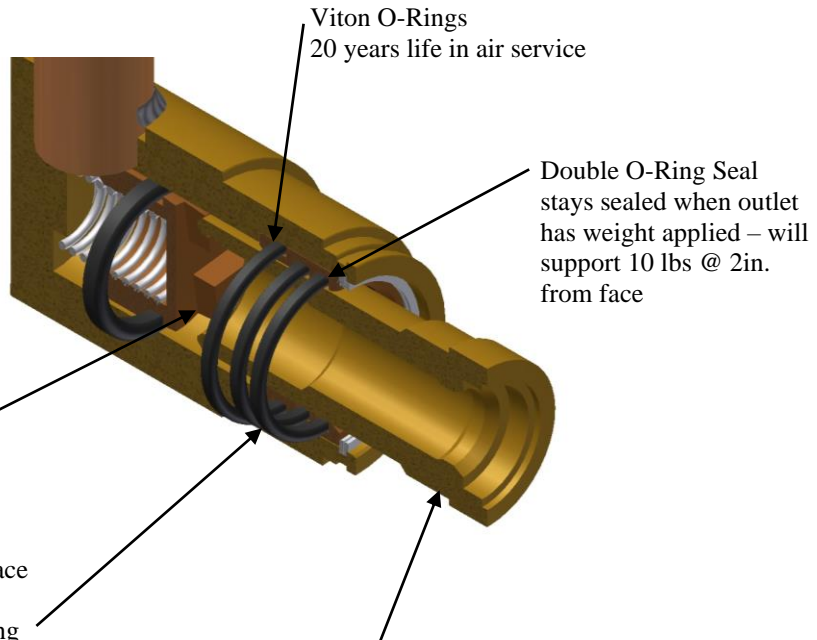


Viton O-Rings  
20 years life in air service

Double O-Ring Seal stays sealed when outlet has weight applied – will support 10 lbs @ 2in. from face

Close tolerance fit supports front body with metal to metal contact to back body

**Typical Competitors Outlet**

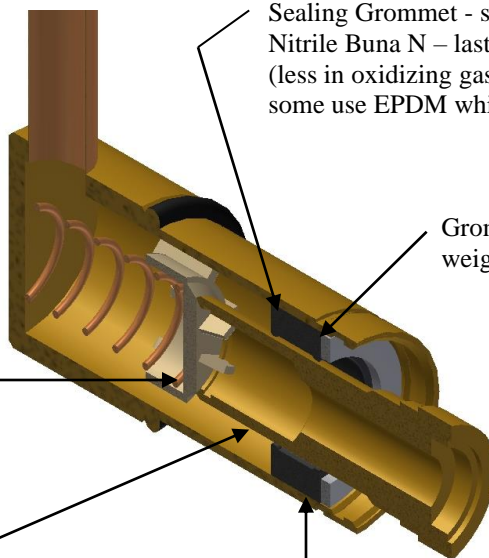

Sealing Grommet - some competitors use a Nitrile Buna N – lasts 3-5 years in air service (less in oxidizing gas services; Oxygen & N<sub>2</sub>O) some use EPDM which lasts somewhat longer

Grommet deforms and leaks when weight is applied to outlet

Plastic Poppet  
Plastic surface seals against rubber grommet

Front body is NOT supported with metal to metal contact deforming grommet and causing LEAKS

To service grommet – gas must be shut off and all parts removed